



αPlastic™
Labware ELISABETH PHARMACON

PRODUCT CATALOG


LABWARE

PCR Tubes
PCR/qPCR Strips
PCR/qPCR Well Plates
Well Plates Seals
Pipette Tips

Content

LABWARE

	Page
PCR Tubes	6
PCR/qPCR Strips	7-9
Tube Opener	10
PCR/qPCR Well Plates	11-13
Well Plates Seals	14
Pipette Tips	15-18

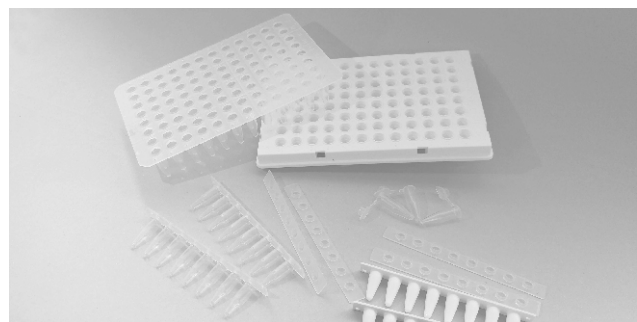



Labware

α Plastic™
Labware ELISABETH PHARMACON

PCR Tubes

PCR 8-Strips

qPCR 8-Strips

PCR 96Well Plates

qPCR 96 Well Plates

Well Plates Seals

Microtubes

Pipette Tips

Tube Opener



αPlastic™ PCR Single Tubes

Features:

- Thin Walled tube for best heat transfer
- Made of USP VI Polypropylene
- Free of DNase, RNase, DNA, PCR Inhibitor, ATP, Endotoxin
- All products are ready to use in PCR
- Compatible with Rotor-Gene and many other cyclers
- Autoclavable
- Evaporation rate below 5%, no mineral oil needed
- Fits almost all thermal cyclers
- Small package can avoid contamination by different researcher
- All items contain tubes with caps

Ref	Description	Features / Functions	Package
αPlastic™ PCR Tube C 0.2 Frosted Flat Cap APTC02FFC-500	0.2 ml PCR Tubes with Flat Cap	Easy labelling on frosted flat cap	500 pcs
αPlastic™ PCR Tube C 0.2 Frosted Flat Cap APTC02FFC-3000	0.2 ml PCR Tubes with Flat Cap	Easy labelling on frosted flat cap	3000 pcs

Parameter	ROX-dependent	No ROX
Machine Brands	ABI, Stratagene Biomtra	Bio-Rad Roche, Eppendorf
Does mastermix include ROX	YES	NO
Ct value determined by	Rn (Reporter RFU/ Rox RFU)	Reporter RFU
Ct value difference among White, Frosted, Clear tubes	White = Frosted = Clear	White > Frosted > Clear

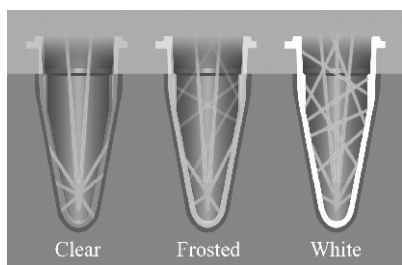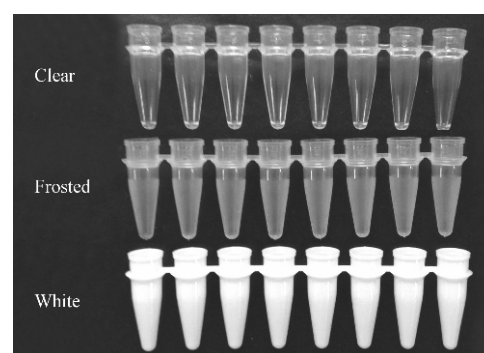

In **no ROX machine**, higher Reporter RFU will determine lower Ct value (the lower Ct the better). The comparison is **White > Frosted > Clear**

In **ROX-dependent machine**, higher Reporter and Passive RFU will get almost the same Rn and determine the same Ct value.
The comparison is **White = Frosted = Clear**

Don't put different type of tube in the same qPCR experiment, that will cause different Ct value, even with the same sample.

COMPATIBILITY CHART (qPCR 8-Strip Tubes)

	AP0801FOC	AP0801WOC	AP0801COC	AP0802COC	AP0802WOC	Note / Holder NO.
Standard Well (0.2 ml, height 21mm.)				✱	✱	
Low Profile Well (0.1 ml, height 16mm.)	✱	✱	✱			
Applied Biosystems (Life Technology)						
7000				✱	✱	
7300				✱	✱	4367033
7500				✱	✱	4367033
7500Fast	✱	✱	✱			4403809
StepOne	✱	✱	✱			4375282
StepOne Plus	✱	✱	✱			4379983
ViiA7 Fast Model (0.1 ml)	✱	✱	✱			4379983
ViiA7 (0.2 ml)				✱	✱	4381850
QuantStudio Fast Model (0.1 ml)	✱	✱	✱			4379983
QuantStudio (0.2 ml)				✱	✱	4381850
QuantStudio 3/5 Fast (0.1 ml)	✱	✱	✱			4379983
QuantStudio 3/5 (0.2 ml)				✱	✱	4381850
3500(xL)	✱	✱	✱			4410232 (Retainer)
(0.1 ml Fast 96 Well)						4410701 (Septa)
3500(xL)				✱	✱	4410230 (Retainer)
(0.2 ml Standard 96 Well)						4410701 (Septa)
Roche						
LightCycler 96	✱	✱	✱			
LightCycler Nano			✱			Bottom Detection
Bio-Rad						
iCycler iQ				✱	✱	
iQ5				✱	✱	
MyiQ				✱	✱	
Chromo4				✱	✱	
Opticon I & II	✱	✱	✱			
mini Opticon	✱	✱	✱			
CFX96	✱	✱	✱			
Stratagene/Eppendorf/Biometra/Sacace/Analytik Jena						
Mx4000/3000P/3005P				✱	✱	
AriaMx	✱	✱	✱			
MyGo PRO			✱			Bottom Detection
Open QPCR			✱			Bottom Detection
MasterCycle ep Realplex				✱	✱	
PikoReal 24 Real-Time PCR	✱	✱	✱			
Toptical Thermocycler				✱	✱	
SaCycler-96				✱	✱	
qTOWER/2.0/2.2				✱	✱	



αPlastic™ PCR 8-Strip Tubes

ENDPOINT

Features:

- Thin Walled tube for best heat transfer
- Made of USP VI Polypropylene
- Free of DNase, RNase, DNA, PCR Inhibitor, ATP, Endotoxin
- All products are ready to use in PCR
- Autoclavable
- Evaporation rate below 5%, no mineral oil needed
- Fits almost all thermal cyclers
- Small package can avoid contamination by different researcher
- All items contain tubes with caps

Ref	Description	Features / Functions	Package
αPlastic™ PCR 8-well strip tube 0.2, Attached Flat Frosted Caps			
AP0802AFFC-120	0.2 ml PCR 8-Strip Tube with Attached Flat Cap	Easy labelling on frosted flat cap Attached cap helps avoid cross-contamination between different tubes Good for storage (10 Bags per 12 strip Tubes)	120 pcs
αPlastic™ PCR 8-well strip tube 0.2, Flat Frosted Caps			
AP0802FFC-120	0.2 ml PCR 8-Strip Tube & Separate Flat Cap	Easy labelling on frosted flat cap Traditional separate cap decrease operating time (10 Bags per 12 strip Tubes)	120 pcs
αPlastic™ PCR 8-well strip tube 0.2, Domed Caps			
AP0802DC-120	0.2 ml PCR 8-Strip Tube & Separate Domed Cap	Dome cap could fit domed hot lid thermal cycler Traditional separate cap decrease operation time (10 Bags per 12 strip Tubes)	120 pcs

αPlastic™ qPCR 8-Strip Tubes

Features:

- Thin Walled tube for best heat transfer
- Made of USP VI Polypropylene
- Free of DNase, RNase, DNA, PCR Inhibitor, ATP, Endotoxin
- All products are ready to use in qPCR
- Evaporation rate below 5%, cost savings on reagents
- Two kind of 8-Strip could fit almost all qPCR machines
- Small package can avoid contamination by different researcher
- Ultra clear cap shows consistent light transmitting rate provide stable data than other brand
- All items contain tubes and caps




Ref	Description	Features / Functions	Package
αPlastic™ qPCR 8-well strip tube LP 0.1F, Optical Caps			
AP0801FOC-120	0.1 ml qPCR 8-Strip Tubes (With Optical Caps)	Low profile Tube height 16 mm Frosted tubes, (10 Bags per 12 strip Tubes)	120 pcs
αPlastic™ qPCR 8-well strip tube LP 0.1W, Optical Caps			
AP0801WOC-120	0.1 ml qPCR 8-Strip White Tubes (With Optical Caps)	Low profile Tube height 16 mm White tubes (10 Bags per 12 strip Tubes)	120 pcs
αPlastic™ qPCR 8-well strip tube LP 0.1C, Optical Caps			
AP0801COC-120	0.1 ml qPCR 8-Strip Clear Tubes (With Optical Caps)	Low profile Tube height 16 mm Ultra clear tubes designed for bottom detection machine like Roche LC Nano, Mygo Pro, Open qPCR (10 Bags per 12 strip Tubes)	120 pcs
αPlastic™ qPCR 8-well strip tube SP 0.2C, Optical Caps			
AP0802COC-120	0.2 ml qPCR 8-Strip Tubes (With Optical Caps)	Standard profile Tube height 21 mm Clear tubes (10 Bags per 12 strip Tubes)	120 pcs
αPlastic™ qPCR 8-well strip tube SP 0.2 W, Optical Caps			
AP0802WOC-120	0.2 ml qPCR 8-Strip White Tubes (With Optical Caps)	Standard profile Tube height 21 mm White tubes (10 Bags per 12 strip Tubes)	120 pcs



αPlastic™ Tube Opener

Features:


- Made of USP VI Polystyrene
- Operating temperature -30 ~ 70 °C
- Multiple functions design for decapping, capping, film scrapping
- Multiple decapping functions could be use in attached 8-strip tubes, single tube
- Transparent design helps observation and operation
- A smooth surface on film scrapping side, avoids causing damage on any adhesive film
- Non-autoclavable, using ethanol spray for sterilization

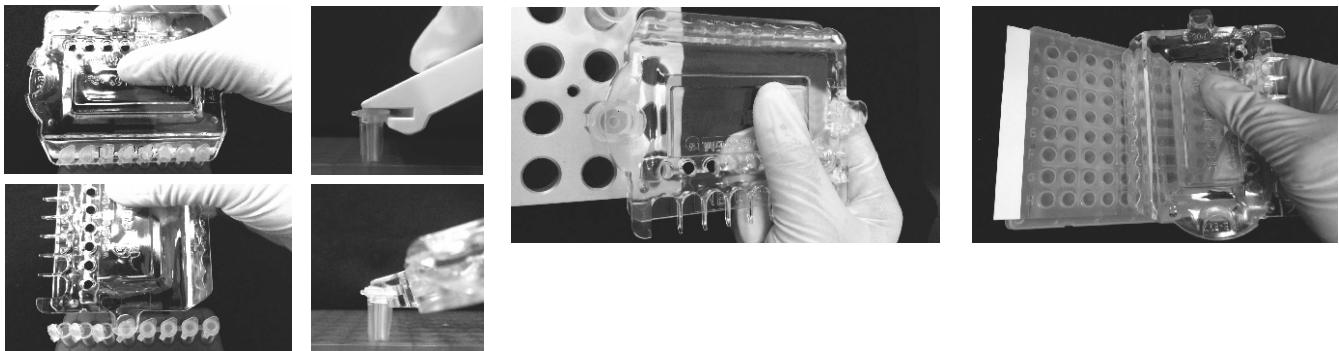

Ref	Package	Compatibility
αPlastic™ Tube Opener		
APCAPPING-1	1 pc	Decapping functions for: αPlastic™ PCR Tube C 0.2 Frosted Flat Cap αPlastic™ PCR 8-well strip tube 0.2, Attached Flat Frosted Caps αPlastic™ PCR 8-well strip tube 0.2, Flat Frosted Caps Capping function for αPlastic™ PCR 8-well strip tube 0.2, Flat Frosted Caps αPlastic™ qPCR 8-well strip tube SP 0.2C, Optical Caps Film Scrapping for αPlastic™ PCR 96 Well-Plates Adhesive Membrane αPlastic™ qPCR 96 Well-Plates Ultra Clear Pressure Membrane αPlastic™ qPCR 96 Well-Plates Ultra Clear Adhesive Membrane αPlastic™ 96 Well-Plates Aluminium Adhesive Membrane
APCAPPING-15	15 pcs	

COMPATIBILITY CHART (PCR & qPCR 96 Well Plates)

	AP9602CNS	AP9602CSSS	AP9602CSS	AP9601FSSS	AP9601FSS	AP9601WSS	Note / Holder NO.
Standard Well (0.2 ml, height 21mm.)	✱	✱	✱				
Low Profile Well (0.1 ml, height 16mm.)				✱	✱	✱	
Applied Biosystems (Life Technology)							
7000		✱					
7300		✱					4367032
7500		✱					4367032
7500Fast				✱			4359652
7900HT Fast Model(0.1 ml)				✱			
7900HT Standard (0.2 ml)		✱					
StepOne Plus				✱			
ViiA7 Fast Model (0.1 ml)				✱			4459846
ViiA7 (0.2 ml)		✱					4459845
QuantStudio Fast Model (0.1 ml)				✱			4459846
QuantStudio (0.2 ml)		✱					4459845
QuantStudio 3/5 Fast (0.1 ml)				✱			
QuantStudio 3/5 (0.2 ml)		✱					
5700/7700		✱	✱				Sequencer
3100/3130 (xL)				✱			4367470 (Retainer)
(0.1 ml Fast 96 Well)							4315933 (Septa)
310/3100/3130 (xL)		✱	✱				4317241 (Retainer)
(0.2 ml Standard 96 Well)							4315933 (Septa)
3700/3730 (xL)				✱			4367469 (Retainer)
(0.1 ml Fast 96 Well)							4315933 (Septa)
3700/3730 (xL)		✱	✱				4410228 (Retainer)
(0.2 ml Standard 96 Well)							4315933 (Septa)
3500(xL)				✱			4409530 (Retainer)
(0.1 ml Fast 96 Well)							4412614 (Septa)
3500(xL)		✱	✱				4410228 (Retainer)
(0.2 ml Standard 96 Well)							4412614 (Septa)
Roche							
LightCycler 480 & 480II						✱	
LightCycler 96						✱	
Bio-Rad							
iCycler iQ	✱		✱				
iQ5	✱		✱				
MyiQ	✱		✱				
Chromo4	✱		✱				
Opticon I & II					✱	✱	
CFX96					✱	✱	
Stratagene/Eppendorf/Biomtra/Sacace/Analytik Jena							
Mx4000/3000P/3005P	✱		✱				
AriaMx					✱	✱	
MasterCycle ep Realplex	✱		✱				
Toptical Thermocycler	✱		✱				
SaCycler-96	✱		✱				
qTOWER/2.0/2.2	✱		✱				

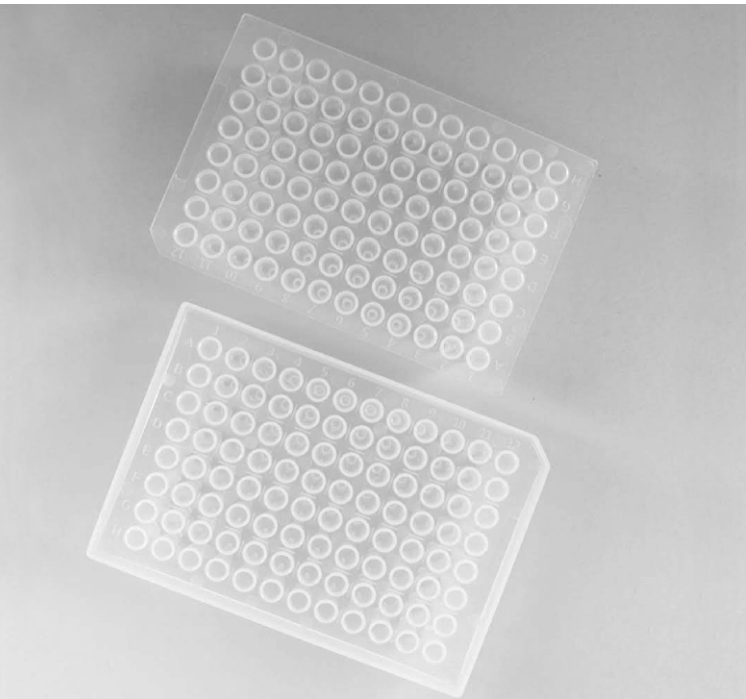


αPlastic™ PCR 96 Well Plates

Features:


- Thin Walled tube for best heat transfer
- Made of USP VI Polypropylene
- Free of DNase, RNase, DNA, PCR Inhibitor, ATP, Endotoxin.
- Evaporation rate below 5%, saving cost on reagents
- Small package can avoid contamination by different researcher



Non-Skirted

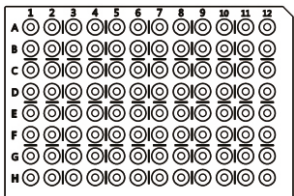

Semi Skirted

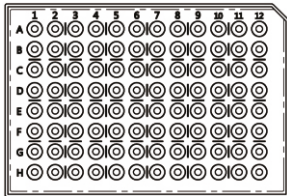


Ref	Package		Description	Features / Functions
αPlastic™ PCR 96 Well-Plates 0.2C Non Skirted				
AP9602CNS-10	10 pcs (1 bag)		0.2 ml PCR 96-Well Plates Non-Skirted	Standard profile Tube height 21 mm Clear tubes, (1 box = 5 bags per 10 plates)
AP9602CNS-50	50 pcs (1 box)			
αPlastic™ qPCR 96 Well-Plates SP 0.2C, Semi Skirted				
AP9602CSS-10	10 pcs (1 bag)		0.2 ml PCR 96 Well Plate Semi Skirted	Standard profile Tube height 21 mm Clear tubes (1 box = 5 bags per 10 plates)
AP9602CSS-50	50 pcs (1 box)			


αPlastic™ qPCR 96 Well Plates

Features:


- Thin Walled tube for best heat transfer
- Made of USP VI Polypropylene
- Free of DNase, RNase, DNA, PCR Inhibitor, ATP, Endotoxin. All products are ready to use in qPCR
- Evaporation rate below 5%, saving cost on reagents
- Low and standard profile with 2 skirted format for almost all qPCR and sequencer machine
- Small package can avoid contamination by different researcher

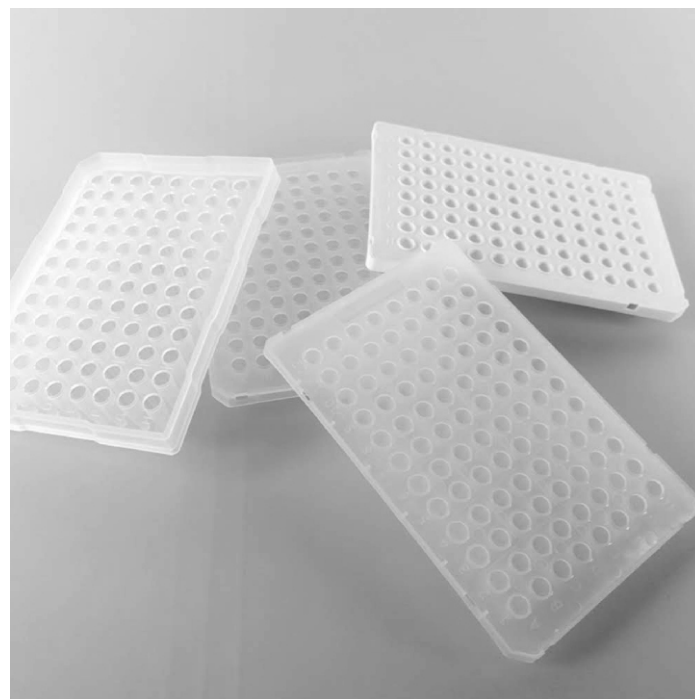

Sub-Semi Skirted

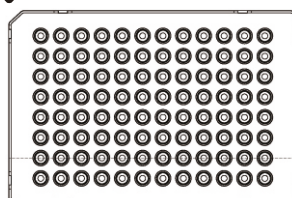

Semi Skirted

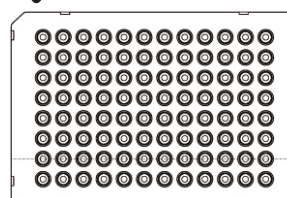


Ref	Package		Description	Features / Functions
αPlastic™ qPCR 96 Well-Plates SP 0.2C, Sub-Semi Skirted				
AP9602CSSS-10	10 pcs (1 bag)		0.2 ml qPCR	Standard profile
AP9602CSSS-50	50 pcs (1 box)		96 Well Plates	Tube height 21 mm
			Sub-Semi Skirted	Clear tubes,
αPlastic™ qPCR 96 Well-Plates LP 0.1F, Sub-Semi Skirted				
AP9601FSSS-10	10 pcs		0.1 ml qPCR	Low profile
AP9601FSSS-50	50 pcs		96 Well Plates	Tube height 16 mm
			Sub-Semi Skirted	Frosted tubes
αPlastic™ qPCR 96 Well-Plates LP 0.1F, Semi Skirted				
AP9601FSS-10	10 pcs		0.1 ml qPCR	Low profile
AP9601FSS-50	50 pcs		96 Well Plates	Tube height 16 mm
			Semi Skirted	Frosted tubes,
αPlastic™ qPCR 96 Well-Plates LP 0.1W Semi Skirted				
AP9601WSS-10	10 pcs		0.1 ml qPCR	Low profile
AP9601WSS-50	50 pcs		96 Well Plates	Tube height 16 mm
			Semi Skirted	White tubes

αPlastic™ Seals


Features:

- Free of detectable DNase, RNase, DNA, PCR Inhibitor, ATP, Endotoxin.
- Products are ready to use in PCR, qPCR or ELISA
- Small package can avoid contamination by different researcher

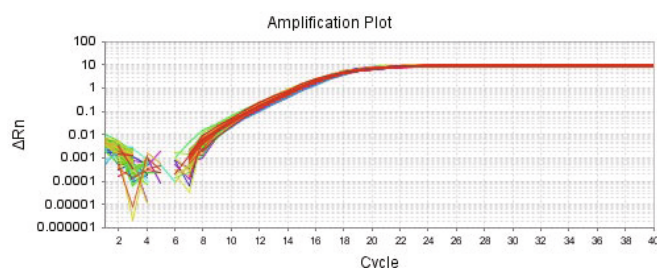

αPlastic™ AP96PMEM-10

The same 96 wells mechanism repeated on ABI StepOnePlus and the data show stable overlap of 96 lines in AP96PMEM-10.

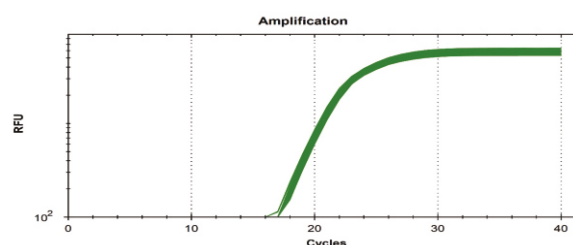

αPlastic™ AP96AMEM-10

The same 96 wells mechanism repeated on CFX-96 and the data show stable overlap of 96 lines in AP96AMEM-10.

Ref	Package	Description	Features / Functions
αPlastic™ PCR 96 Well-Plates Adhesive Membrane			
AP96MEM-10	10 pcs (1 bag)	PCR 96 Well Plate	Operating Temperature -40°C ~ 120°C PCR use Strong adhesive
AP96MEM-100	100 pcs (1 box)	Sealing Membrane	
αPlastic™ qPCR 96 Well-Plates Ultra Clear Pressure Membrane			
AP96PMEM-10	10 pcs (1 bag)	qPCR 96 Well Plate	Operating Temperature -40°C ~ 120°C Real-Time PCR(qPCR) use Strong pressure sensitive adhesive and high light transmitting rate provides more table data
AP96PMEM-100	100 pcs (1 box)	Optical Sealing Membrane (pressure-sensitive adhesive film) equivalent: ABI 4311971	
αPlastic™ qPCR 96 Well-Plates Ultra Clear Adhesive Membrane			
AP96AMEM-10	10 pcs (1 bag)	qPCR 96 Well Plate	Operating Temperature -40°C ~ 120°C Real-Time PCR(qPCR) use Strong adhesive Ultra clear membrane shows consistent and high light transmitting rate, provides more stable data
AP96AMEM-100	100 pcs (1 box)	Optical Sealing Membrane (adhesive film) equivalent: BIO-RAD MSB-1001	
αPlastic™ 96 Well-Plates Aluminium Adhesive Membrane			
AP96ALMEM-10	10 pcs (1 bag)	Aluminium Membrane	Operating Temperature -80°C ~ 120°C PCR or Storage use Strong adhesive MB-ASM is pierceable
AP96ALMEM-100	100 pcs (1 box)		
αPlastic™ 96 Well-Plates ELISA CLEAR Adhesive Membrane			
AP96EMEM-100	100 pcs (1 box)	Aluminium Membrane	Clear Storage membrane Elisa/EIA or various storage use Operating Temperature -40°C ~ 120°C

αPlastic™ Pipette Tips

We are offering a wide range of reliable **αPlastic Pipette Tips** for confident use in your laboratory experiments.


Our **Low Retention tips** are produced with technology that creates an ultrahydrophobic non-wettable surface that limits the adhesion of cells, DNA/RNA, proteins, and other viscous materials. The production process does not use surface coatings that can lead to sample contamination or PCR inhibition. Low retention properties ensure full volume retrieval of your sample. The pipetting accuracy and handling precision is up to 10 times better with viscous samples than

other standard polypropylene tips. As a result of this our tips offer excellent sample recovery while increasing the accuracy and reproducibility of microvolume pipetting in your experiments.

αPlastic Pipette Tips have a clear graduation marking on the outside of the tip to help visually confirm pipette settings. These sterile tips are also produced with superior filter barrier that protects your pipette and samples from aerosol cross-contamination during PCR. Universally-fitting tips are designed to fit practically all single and multichannel pipettes. This compatibility ensures that you can use them with most leading pipette brands (see the Compatibility Table).

Finally all of pipette tips listed in this catalog are RNase-, DNase- and Endotoxin free and produced with high quality, non-cytotoxic medical grade, 100% radiation stabilized virgin polypropylene. These resins are pre-tested from molded samples to determine whether any leachable compounds that could interfere with your experiments are present.


Our highest priority is to offer you reliable tips with excellent quality and for competitive prices. This makes our **αPlastic Pipette Tips** your No.1 choice for your experiments.

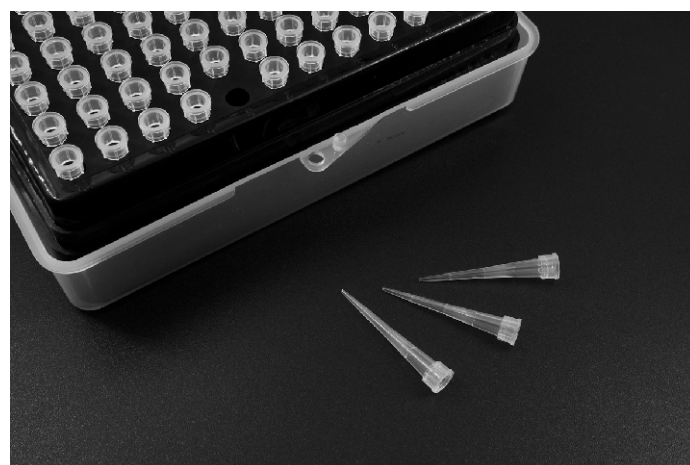


αPlastic™ 10µl Universal tip, Sterile, Low retention, Filter

Universal **αPlastic 10 µl Pipette Tips** for a wide range of biomolecular applications

Ref	Package	Description
NT-10XN-S-L-F-R	96 tips (min. 10x96)	Extended Length
C-10-S-L-F-R	96 tips (min. 10x96)	Standard

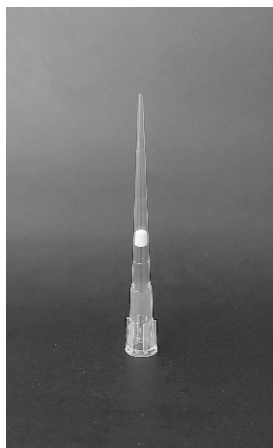

NT-10XN-S-L-F-R

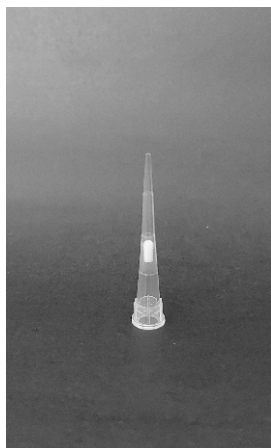


C-10-S-L-F-R

Key Features


- High Quality polypropylene (PP)
- Super hydrophobic internal surface for maximum recovery of aspirate
- Standard size and Extra Long and Narrow tip to reach the bottom of 1.5 ml micro centrifuge tubes
- Ideal for both single-channel and multi-channel pipettes
- High-density filter to prevent cross-contamination
- Aerosol resistant
- Universal compatibility
- Low retention for precise pipetting
- Free of detectable RNase, DNase, DNA, and pyrogens
- Graduation marks for accurate volumes transfer
- Sterile packaging

αPlastic™ 200µl Universal tip, Sterile, Low retention, Filter

Universal αPlastic 200 µl Pipette Tips for a wide range of biomolecular applications

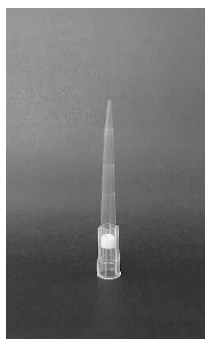

NT-200-S-L-F-R

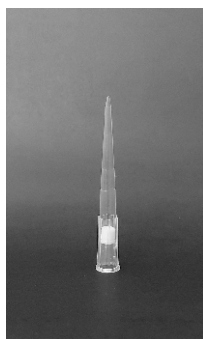

C-200-S-L-F-R


Ref	Package	Description
NT-200-S-L-F-R	96 tips (min. 10x96)	Standard
C-200-S-L-F-R	96 tips (min. 10x96)	Standard

Key Features

- High Quality polypropylene (PP)
- Super hydrophobic internal surface for maximum recovery of aspirate
- Ideal for both single-channel and multi-channel pipettes
- High-density filters to prevent cross-contamination
- Aerosol resistant
- Universal compatibility
- Low retention for precise pipetting
- Free of detectable RNase, DNase, DNA, and pyrogens
- Graduation marks for accurate volumes transfer
- Sterile packaging

αPlastic™ 1000µl Universal tip, Sterile, Low retention, Filter

Universal αPlastic 1000 µl Pipette Tips for a wide range of biomolecular applications

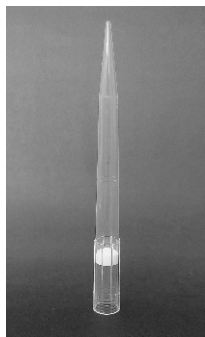

NT-1000-S-L-F-R


Ref	Package	Description
NT-1000-S-L-F-R	96 tips (min. 6x96)	Standard

Key Features

- High Quality polypropylene (PP)
- Super hydrophobic internal surface for maximum recovery of aspirate
- Ideal for both single-channel and multi-channel pipettes
- High-density filters to prevent cross-contamination
- Aerosol resistant
- Universal compatibility
- Low retention for precise pipetting
- Free of detectable RNase, DNase, DNA, and pyrogens
- Graduation marks for accurate volumes transfer
- Sterile packaging

αPlastic™ 1250µl Universal tip, Sterile, Low retention, Filter

Universal αPlastic 1250 µl Pipette Tips for a wide range of biomolecular applications

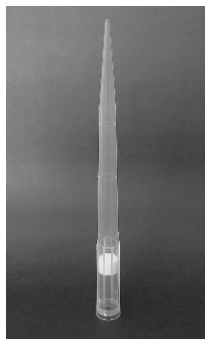

C-1250-S-L-F-R

Ref	Package	Description
C-1250-S-L-F-R	96 tips (min. 8x96)	Extended Length

Key Features

- High Quality polypropylene (PP)
- Super hydrophobic internal surface for maximum recovery of aspirate
- Extra Long sized tip to reach the very bottom of each tube
- Increased space between tip and filter for fast repetitive pipetting of large volumes
- Ideal for single-channel pipettes only
- High-density filters to prevent cross-contamination
- Aerosol resistant
- Universal compatibility
- Low retention for precise pipetting
- Free of detectable RNase, DNase, DNA, and pyrogens
- Graduation marks for accurate volumes transfer
- Sterile packaging

COMPATIBILITY CHART (αPlastic™ Pipette Tips)

Certified DNase, RNase, Pyrogen, Endotoxin, Human DNA & PCR Inhibitor Free 		NT-10XN-S-L-F-R	NT-200-S-L-F-R	NT-1000-S-L-F-R
		10 µl Filter tips	200 µl Filter tips	1000 µl Filter tips
BIOHIT	mLine Advanced Digital Pipette			
	BM-200 Biohit mLine manual pipettor, 1 Ch. 20-200ul		●	
	BM-1000 Biohit mLine manual pipettor, 1 Ch. 100-1000ul			●
	Proline® Plus Single Channel Pipettes			
	BU-200 Proline Plus 1 channel pipettor 20-200ul		●	
	BU-1000 Proline Plus 1 channel pipettor 100-1000ul			●
	Proline Mechanical Pipette			
	BH-200 Biohit Proline Manual Pipettor, 1 ch., 20-200ul		●	
	BH-1000 Biohit Proline Manual Pipettor, 1 ch., 100-1000ul			●
	eLine Advanced Electronic Pipette			
	BL-120A Biohit Eline Electronic Pipettor 1 ch. 5-120ul w/Adapter		●	
	Proline Electronic Pipette			
SARTORIUS	BP-500 Biohit Proline Electronic Pipettor. Single channel, 10-500ul in 5ul increments			●
	BP-1000 Biohit Proline Electronic Pipettor. Single channel, 50-1000ul in 10ul increments			●
	Tacta Mechanical Pipette			
	LH-729010 Tacta Mechanical 0.1-3 µl single channel	●		
	LH-729020 Tacta Mechanical 0.5-10 µl single channel	●		
	LH-729060 Tacta Mechanical 20-200 µl single channel		●	
	LH-729070 Tacta Mechanical 100-1000 µl single channel			●
	Picus			
	735021 Picus 0.2-10 ul single channel	●		
	735081 Picus 50-1000 ul single channel			●
	Picus Nxt			
	LH 735021 Picus 0.2-10 ul single channel	●		
	LH 735081 Picus 50-1000 ul single channel			●
CAPP	CU-200 CappAero 1 ch. pipette, 20-200 ul, Ecoline, Yellow		●	
	CU-1000 CappAero 1 ch. pipette, 100-1000 ul, Ecoline, Blue			●
EPENDORF	Reference Model 2000 PhysioCare Manual Pipette			
	E2000-10 Eppendorf Reference 2000 Model 1 channel, 1-10ul Air-Displacement M.	●		
	E2000-200 Eppendorf Reference 2000 Model 1 channel, 50-200ul Air-Displacement M.		●	
	E2000-1000 Eppendorf Reference 2000 Model 1 channel, 100-1000ul Air-Displacement M.			●
	Research Model 2100 PhysioCare Manual Pipette			
	E2100-2.5 Eppendorf Research 2100 1 channel, 0.1-2.5ul Air-Displacement M.	●		
	E2100-10 Eppendorf Research 2100 1 channel, 1-10ul Air-Displacement M.	●		
	E2100-200 Eppendorf Research 2100 1 channel, 20-200ul Air-Displacement M.		●	
	E2100-1000 Eppendorf Research 2100 1 channel, 100-1000ul Air-Displacement M.			●
	Research Plus Single Channel Pipette			
	Research Plus 0.5-10 ul	●		●
	Research Plus 2-20 ul (yellow)		●	
	Research Plus 10-100 ul		●	
	Research Plus 20-200 ul		●	
	Research Plus 100-1000 ul			●
	Research Plus Multi Channel Pipette			
	Research Plus 0.5-10 ul	●		
	Research Plus 10-100 ul		●	
	Research Plus 30-300 ul		●	
	Research Pro Electronic Pipette			
	EE-1000 Eppendorf Research pro 1 channel, 100-1000ul Electronic Pipette			●
	EE-10	●		
	EE12-10	●		
	Xplorer / Xplorer Plus Single Channel			
	0.5-10 ul	●		
	5-100 ul		●	
	15-300 ul		●	
	50-1000 ul			●
	Xplorer / Xplorer Plus Multi Channel			
	0.5-10 ul	●		
	5-100 ul		●	
	15-300 ul		●	
	50-1200 ul			●
Fisher Scientific	FisherBrand Finnpiptette I Digital Pipette			
	FFI-200 FisherBrand Finnpiptette I 20-200 ul		●	
	FFI-1000 FisherBrand Finnpiptette I 100-1000 ul			●

- Complete compatibility table of all tips can be found at www.elisabeth.cz in the tips section

COMPATIBILITY CHART (αPlastic™ Pipette Tips)


Certified DNase, RNase, Pyrogen, Endotoxin, Human DNA & PCR Inhibitor Free				
		NT-10XN-S-L-F-R	NT-200-S-L-F-R	NT-1000-S-L-F-R
		10 µl Filter tips	200 µl Filter tips	1000 µl Filter tips
GILSON	Pipetman Traditional Shaft Pipette			
	P-10 Gilson Pipetman 1 channel Air-Displacement Micropipette, 1 - 10 ul	•		
	P-200 Gilson Pipetman 1 channel Air-Displacement Micropipette, 50 - 200 ul		•	
	P-1000 Gilson Pipetman 1 channel Air-Displacement Micropipette, 100 - 1000 ul			•
	Pipettman Ultra			
Labnet HTL	U10 1-10 ul	•		
	U1000 200-1000 ul			•
	Labpette Discovery Digital Pipette			
Labnet HTL	LN-2 Labnet Labpette Discovery 0.2 to 2ul autoclavable	•		
	LN-200 Labnet Labpette Discovery 20 to 200ul autoclavable		•	
	LN-1000 Labnet Labpette Discovery 100 to 1000ul autoclavable			•
NICHIRYO	Nichipet EX Single Channel Pipette			
	NPX-2 Nichipet EX Autoclavable/UV stable digital micropipette 0.2-2ul	•		
	NPX-10 Nichipet EX Autoclavable/UV stable digital micropipette 1-10ul	•		
	NPX-1000 Nichipet EX Autoclavable/UV stable digital micropipette 100-1000ul			•
	Nichipet EX 2 Single Channel Pipette			
NICHIRYO	NPX2-2 Nichipet EX Autoclavable/UV stable digital micropipette 0.1-2ul	•		
	NPX2-10 Nichipet EX Autoclavable/UV stable digital micropipette 10-10ul	•		
	NPX2-200 Nichipet EX Autoclavable/UV stable digital micropipette 20-200ul		•	
RAININ*	NPX2-1000 Nichipet EX Autoclavable/UV stable digital micropipette 100-1000ul			•
	Electronic Pipette			
	ED1-200R Rainin EDP 1 channel Electronic Pipette, 20 - 200 ul		•	
RAININ*	ED2-1000R Rainin EDP 2 1 channel Electronic Pipette, 100 - 1000 ul			•
	ED-1000R Rainin EDP Plus channel Electronic Pipette, 100 - 1000 ul			•
	Acura 825 Single Channel Digital Pipette			
SOCOREX	S825-200 Socorex Acura® 825 Single-Channel Variable Volume Pipettor 20 - 200		•	
	S825-1000 Socorex Acura® 825 Single-Channel Variable Volume Pipettor 100 - 1000			•
	Transferpette			
BRAND	704744 Transferpette 0.5-10 ul	•		
	704778 Transferpette 20-200 ul		•	
	704780 Transferpette 100-1000 ul			•

*Note : Our Tips are compatible with Rainin Classic Pipettes

αPlastic Tips

Length

C-10-S-L-F-R	31.65 mm
NT-10Nx-S-L-F-R	45.80 mm
C-200-S-L-F-R	50.55 mm
NT-200-S-L-F-R	53.40 mm
NT-1000-S-L-F-R	83.50 mm
C-1250-S-L-F-R	102,01 mm

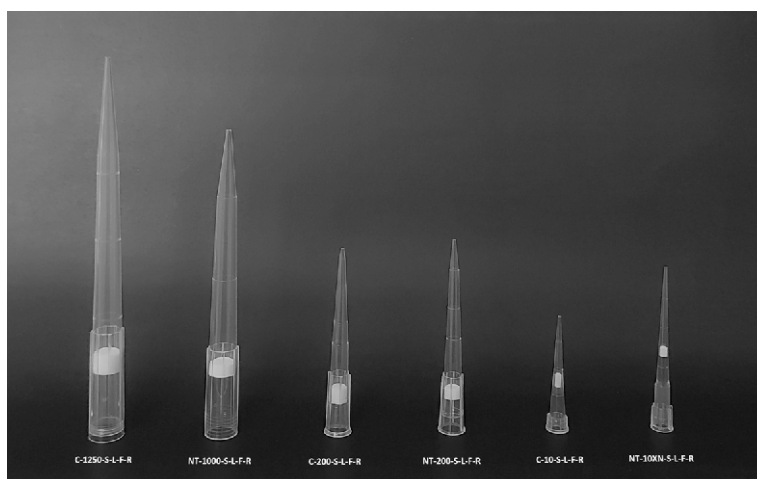

- Complete compatibility table of all tips can be found at www.elisabeth.cz in the tips section

www.elisabeth.cz | www.eligene.com | www.elizyme.com

Manufacturer:

ELISABETH PHARMACON, spol. s r.o.
Rokycanova 4437/5, 615 00 Brno-Zidenice
phone: +420 542 213 851
e-mail: info@elisabeth.cz
www.elisabeth.cz

Distributor:

